


NexSite™

Tunneled Step Tip Hemodialysis Catheter


*Designed to facilitate healing
at the catheter exit site*

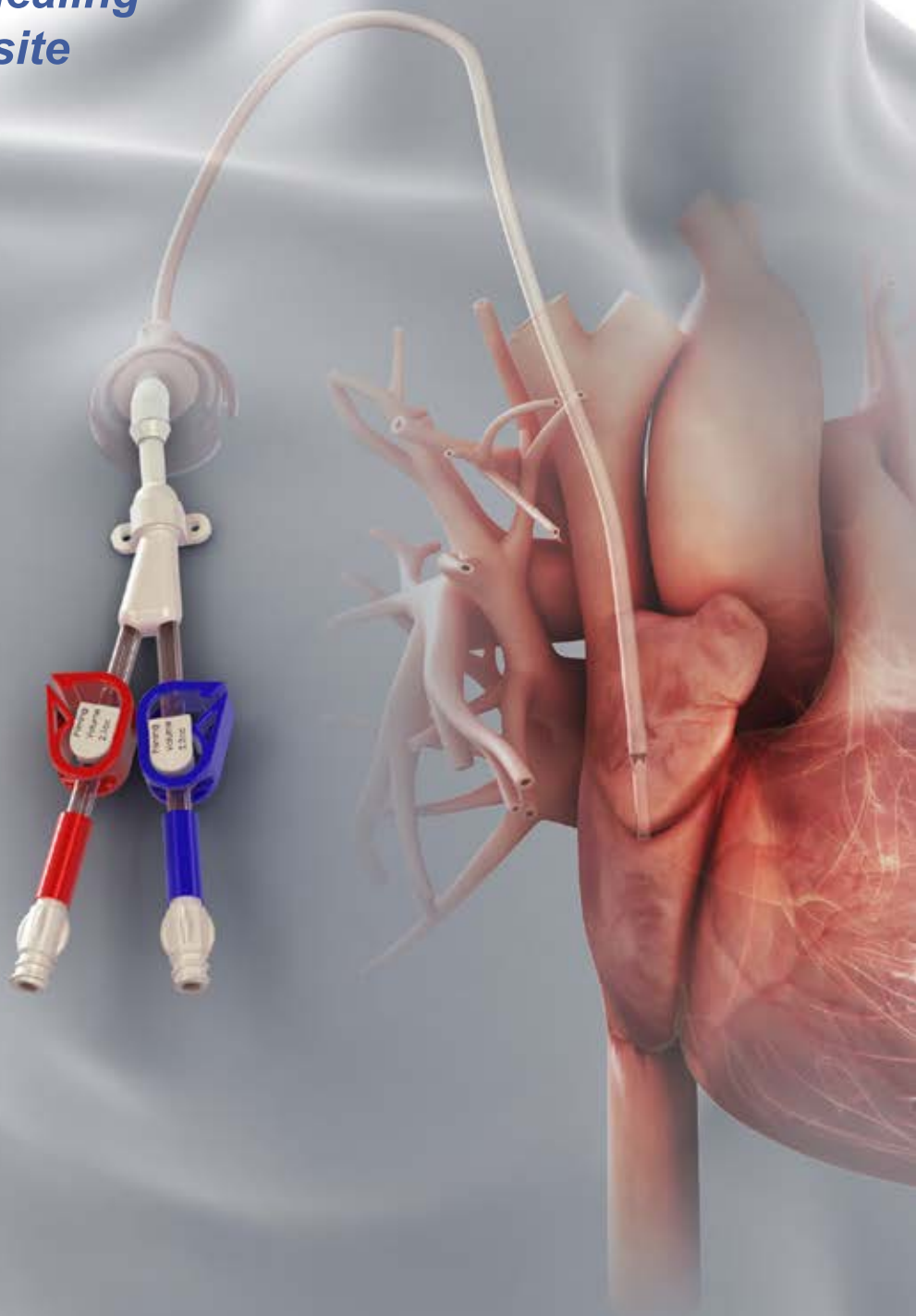

MARVAOMEDICAL

The NexSite™ DISC™ (Dermal Ingrowth Support Collar) features a biomaterial tissue ingrowth scaffold and a guide that assists with the direction of the catheter. The polyester biomaterial scaffolds on the catheter shaft and DISC are aligned to facilitate tissue ingrowth, which is considered important in reducing the source of extraluminal infection in patients requiring long term catheterization.

NexSite™ ESM™ (Exit Site Management) Technology

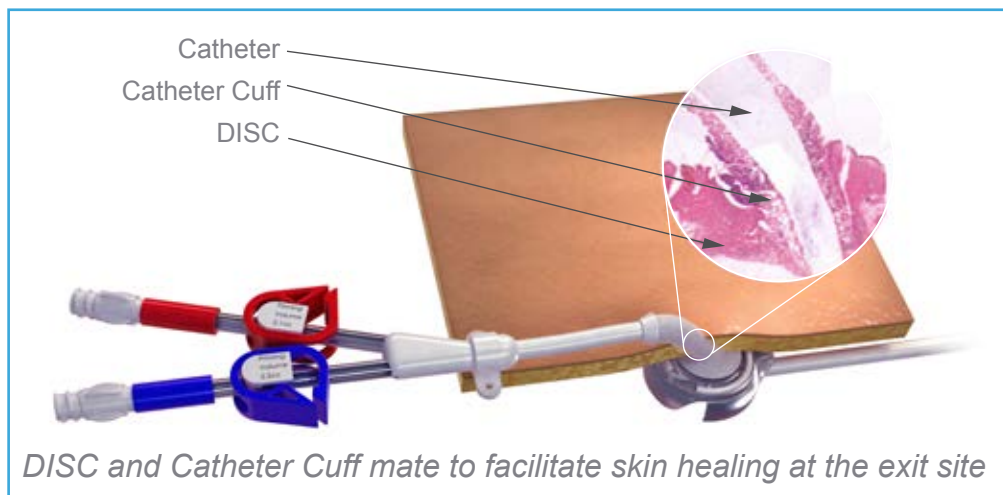

NexSite First in Man Study:

NexSite facilitates healing at the catheter exit site to establish a natural biological seal between the skin and catheter.



NexSite Hemodialysis Catheter Description:

- Radiopaque 15.5Fr step-tip dual-lumen catheter with bifurcation hub, extension tubes and luer connectors
- Manufactured from extruded polyurethane
- Available in 24cm, 28cm, 32cm, 36cm, 40cm and 55cm lengths
- All catheters are available with a transcutaneous polyester cuff and DISC

Package: Kit - 5 each/carton

Product Codes:

24cm NexSite Step Tip Hemodialysis Catheter.....NEXHD1552401
 28cm NexSite Step Tip Hemodialysis Catheter.....NEXHD1552801
 32cm NexSite Step Tip Hemodialysis Catheter.....NEXHD1553201
 36cm NexSite Step Tip Hemodialysis Catheter.....NEXHD1553601
 40cm NexSite Step Tip Hemodialysis Catheter.....NEXHD1554001
 55cm NexSite Step Tip Hemodialysis Catheter.....NEXHD1555501



Marvao Medical Devices, Ltd

liBC at GMIT
 Dublin Road, Galway, Ireland
 p +353 (0)91 759 301
www.marvaomedical.com